



## Canary-Safe™ Exhaust Filters

Exhaust filters are an essential part of the Canary Safety system. Canary-Safe exhaust filters contain an activated carbon, multi-component granulate material that offers optimum filtration for nearly all solvent vapors. With a specific surface area of 1,2000m<sup>2</sup>/g, it absorbs 99% of all vapor. The filter material contains components that prevent sticking or clumping to ensure peak performance, and prevents particulate material from evaporating from your solvents. Filters come with a Timer/Change Indicator at no additional cost.

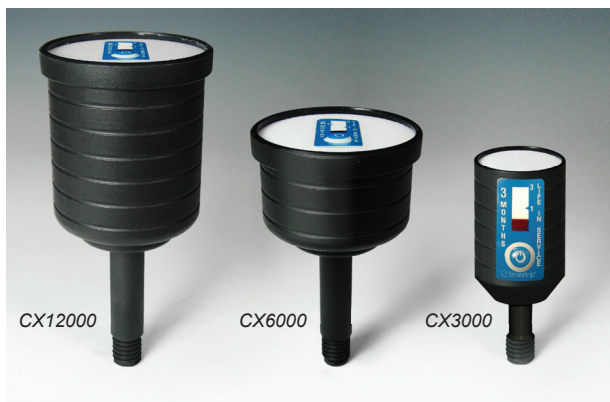

### Timer/Change Indicator

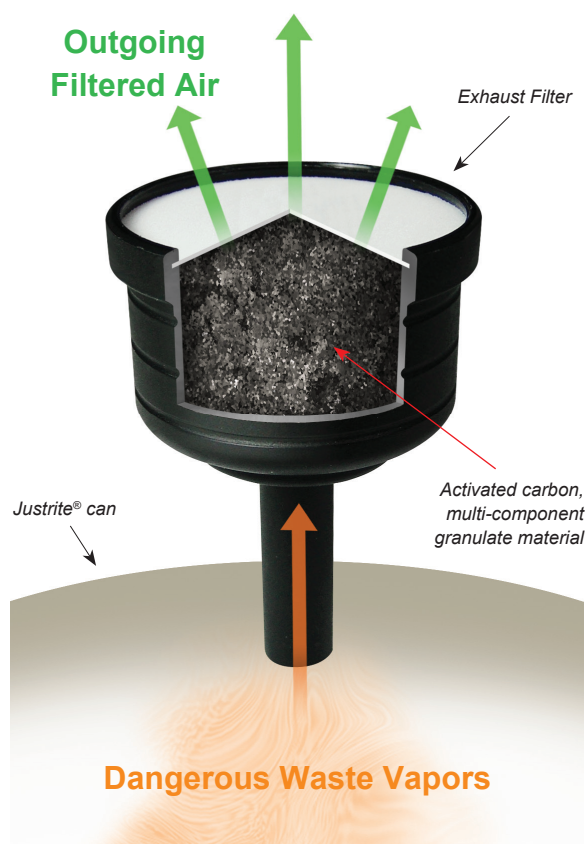

Push button to activate



Replace when service life expires

### Canary-Safe™ Vapor Exhaust Filters

Cat. No.	Size	Description	Replacement
CX3000	S	Exhaust Filter, for containers < 5L	3 Months
CX6000	M	Exhaust Filter, for containers < 20L	6 Months
CX12000	L	Exhaust Filter, for containers > 20L	12 Months
CX6010	M	Exhaust Filter with H <sub>2</sub> SO <sub>4</sub> (sulfuric acid) for Volatile Bases	6 Months
CX6020	M	Exhaust Filter with KOH (potassium hydroxide) for Volatile Acids	6 Months

### Exhaust Filter Adapters for JustRite® Waste Containers

Cat. No.	Description	Qty
CJX01	Exhaust Filter Fitting/ <b>Long</b> Adapter for Stainless Steel Port on JustRite Waste Container	Each
CJX02	Exhaust Filter Fitting/ <b>Short</b> Adapter for Gray PEEK Port on JustRite Waste Container	Each

