


*Ideal for use in
automated liquid handlers*


Storage/Assay 24 & 48 Well Plates



Cat. No.
24108

| Format | Partitions | Qty |
|----------------|-------------------|-----|
| 24-Well | 6 Columns, 4 Rows | 25 |

Deep Well PP Square Well Storage/Reaction Microplate, 10mL/well, 240mL max, Height: 44mm, Well Bottom Shape: **Pyramid**



Cat. No.
24727


| Format | Partitions | Qty |
|----------------|-------------------|-----|
| 24-Well | 6 Columns, 4 Rows | 25 |

Deep Well PP Square Well Storage/Reaction Microplate, 10mL/well, 240mL max, Height: 44mm, Well Bottom Shape: **Round**



Cat. No.
24777


| Description | Qty |
|---|-----|
| Silicone/PTFE Cap Mat for 24-well Collection Plates | 5 |



Cat. No.
48685

| Format | Partitions | Qty |
|----------------|-------------------|-----|
| 48-Well | 6 Columns, 8 Rows | 25 |


Deep Well PP Rectangle Well Collection Plate, 5mL/well, Height: 44mm, Well Bottom Shape: **Pyramid**

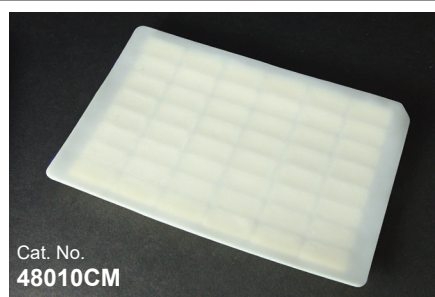

Cat. No.
48010CM

| Description | Qty |
|---|-----|
| Pierceable Santoprene Cap Mat for 48-Well Plates, White | 100 |


Reagent Reservoirs

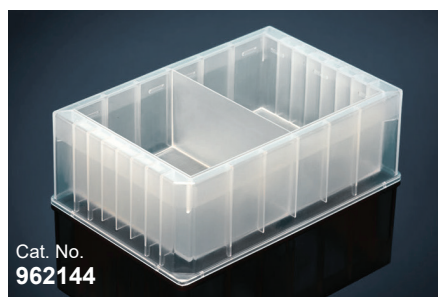

Cat. No.
962144

| Format | Partitions | Qty |
|----------------|------------|-----|
| Columns | 2 Columns | 25 |


Deep Well PP Reservoir, 144mL/Column, 288mL max, Height: 44mm, Well Bottom Shape: **Pyramid**



Cat. No.
960395

| Format | Partitions | Qty |
|----------------|------------|-----|
| Columns | 3 Columns | 25 |

Deep Well PP Reservoir, 95mL/Column, 285mL max, Height: 44mm, Well Bottom Shape: **Pyramid**




Cat. No.
96733

| Format | Partitions | Qty |
|----------------|------------|-----|
| Columns | 4 Columns | 25 |

Deep Well PP Reservoir, 73mL/column, 292mL max, Height: 44mm, Well Bottom Shape: **Pyramid**

Reagent Reservoirs (cont.)

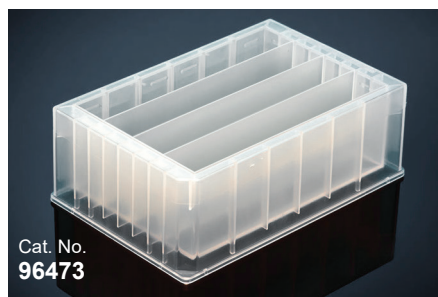


| Format | Partitions | Qty |
|-------------|------------|-----|
| Rows | 4 Rows | 25 |

Deep Well PP Reservoir, 73mL/row, 292mL max, Height: 44mm, Well Bottom Shape: Pyramid




| Format | Partitions | Qty |
|-------------|------------|-----|
| Rows | 8 Rows | 25 |

Low Profile PP Reservoir, 10.6mL/Row, 85mL max, Height: 19mm, Well Bottom Shape: Pyramid

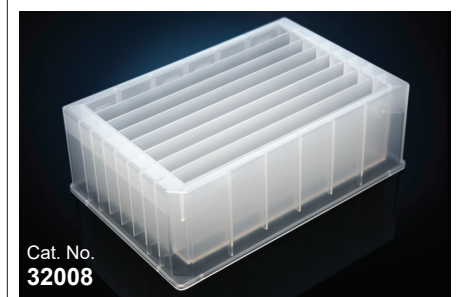

| Format | Partitions | Qty |
|-------------|------------|-----|
| Rows | 8 Rows | 25 |

Deep Well PP Reservoir, 32mL/row, 256mL max, Height: 44mm, Well Bottom Shape: Pyramid

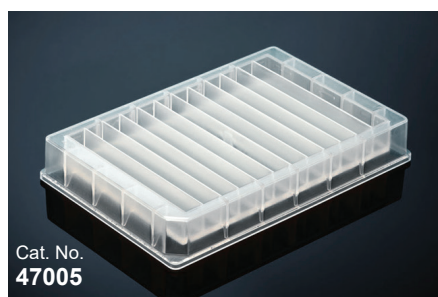


| Format | Partitions | Qty |
|----------------|------------|-----|
| Columns | 12 Columns | 25 |

Low Profile PP Reservoir, 7mL/column, 84mL max, Height: 19mm, Well Bottom Shape: Pyramid




| Format | Partitions | Qty |
|----------------|------------|-----|
| Columns | 12 Columns | 25 |

Deep Well PP Reservoir, 21mL/column, 252mL max, Height: 44mm, Well Bottom Shape: Pyramid

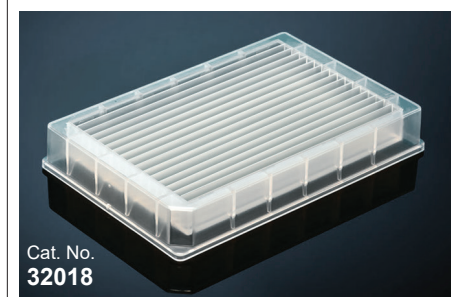

| Format | Partitions | Qty |
|-------------|------------|-----|
| Rows | 16 Rows | 25 |

Low Profile PP Reservoir, 4.9mL/Row, 78mL max, Height: 19mm, Well Bottom Shape: Pyramid

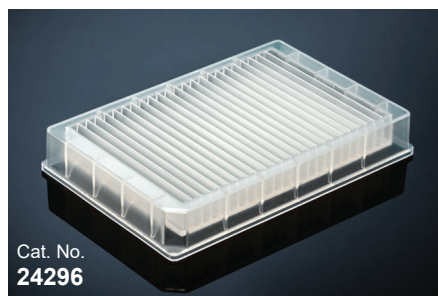


| Format | Partitions | Qty |
|----------------|------------|-----|
| Columns | 24 Columns | 25 |

Low Profile PP Reservoir, 3.3mL/column, 78mL max, Height: 19mm, Well Bottom Shape: Pyramid



| Format | Partitions | Qty |
|----------------|--------------------|-----|
| Special | 12 Columns, 2 Rows | 25 |

Low Profile PP Reservoir, 3.5mL/column, 84mL max, Height: 19mm, Well Bottom Shape: Pyramid

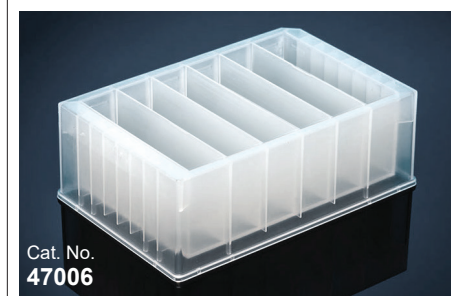


| Format | Partitions | Qty |
|----------------|-------------------------|-----|
| Columns | 6 Columns, 24-Well Base | 25 |


Deep Well PP Reservoir, 47mL/column, 282mL max, Height: 44mm, Well Bottom Shape: Pyramid

Reservoir/Collection Plate Lids

- Universal, 96-Well Robot Compatible
- Works great on standard deep well plates

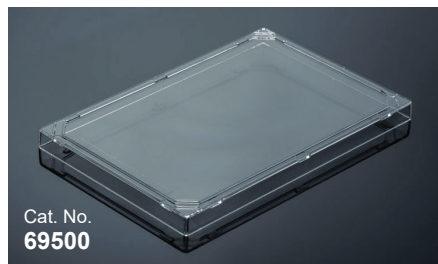


| Cat. No. | Description | Qty |
|----------|-----------------------------|-----|
| 69500 | Reservoir Lid, Clear | 100 |
| 69500BK | Reservoir Lid, Black | 100 |

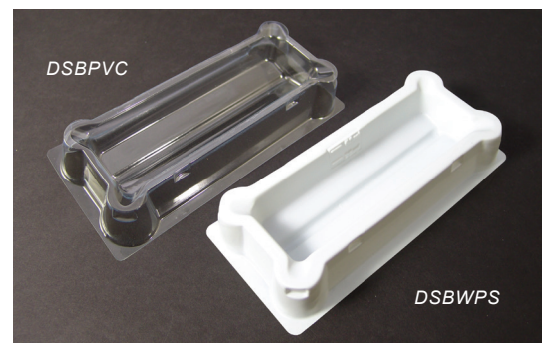

Disposable Basins

| Cat. No. | Description | Bottom | Plate Volume | Qty |
|----------|---|----------|--------------|-----|
| DSBWPS | Disposable Basin, Angled Sides, White PS | V-Bottom | 55mL | 100 |
| DSBPVC | Disposable Basin, Angled Sides, Clear PVC | V-Bottom | 55mL | 100 |

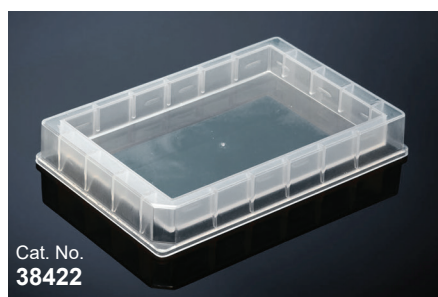

Cat. No.
38422

| Format | Partitions | Qty |
|-------------|------------|-----|
| Open | None | 25 |


Low Profile PP Reservoir, 90mL max Volume, Height: 19mm, Well Bottom Shape: Flat



Cat. No.
47170

| Format | Partitions | Qty |
|-------------|----------------------|-----|
| Open | None, 12 Column Base | 25 |


Half Height PP Reservoir, no partitions, 170mL max, Height: 30.6mm, Column Shape: Pyramid

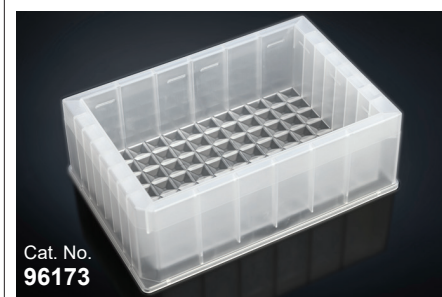

Cat. No.
96173

| Format | Partitions | Qty |
|-------------|--------------------|-----|
| Open | None, 96-Well Base | 25 |


Deep Well PP Reservoir, 287mL max Volume, Height: 44mm, Well Bottom Shape: Pyramid

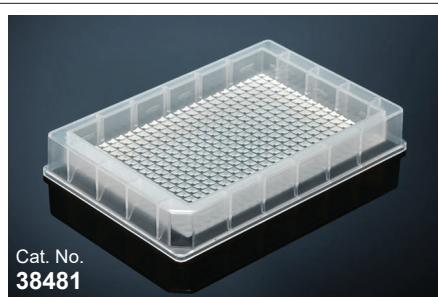

Cat. No.
38481

| Format | Partitions | Qty |
|-------------|---------------------|-----|
| Open | None, 384 Well Base | 25 |


Low Profile PP Reservoir, 92mL max Volume, Height: 19mm, Well Bottom Shape: Pyramid

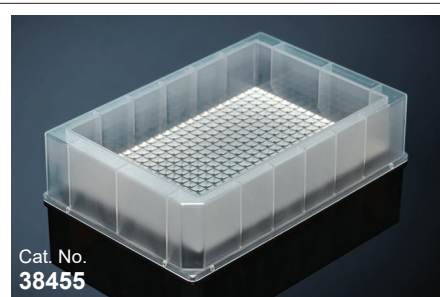

Cat. No.
38455

| Format | Partitions | Qty |
|-------------|---------------------|-----|
| Open | None, 384-Well Base | 25 |


Half Height PP Reservoir, 170mL max Volume, Height: 30.6mm, Well Bottom Shape: Flat

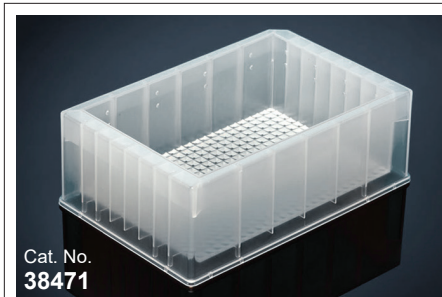

Cat. No.
38471

| Format | Partitions | Qty |
|-------------|---------------------|-----|
| Open | None, 384-Well Base | 25 |

Deep Well PP Reservoir, 282mL max Volume, Height: 44mm, Well Bottom Shape: Pyramid

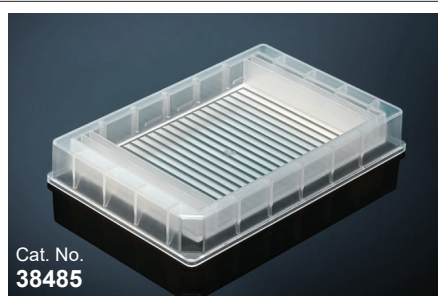

Cat. No.
38485

| Format | Partitions | Qty |
|----------------|----------------------------|-----|
| Special | 4 Controls, 24 Column Base | 25 |

Low Profile PP Reservoir, 3.25mL/control, 86.3mL main, Height: 19mm, Well Bottom Shape: Pyramid



Cat. No.
38444

| Format | Partitions | Qty |
|----------------|---------------------------|-----|
| Special | 4 Controls, 384-Well Base | 25 |

Half Height PP Reservoir, 2.6mL/control, 155mL main, Height: 30.6mm, Well Bottom Shape: Pyramid



Cat. No.
38477

| Format | Partitions | Qty |
|----------------|---------------------------|-----|
| Special | 2 Controls, 384-Well Base | 25 |

Deep Well PP Reservoir, 7mL & 2mL/control, 273mL main, Height: 30.6mm, Bottoms: Pyramid

V-Troughs



96013



96012

NEW!



55001



Cover



| Cat. No. | Description | Qty |
|----------|---|-----|
| 96013 | 12 Column V-Trough, Clear PETG | 5 |
| 96012 | 8 Row V-Trough, White Polystyrene | 5 |
| 55001 | Polypropylene dual reservoir with Lid, 50mL on one side, flip over - 12, 5mL wells on opposite side | 25 |



Vortexer Multipurpose Mixer

This sleek and modern designed vortex mixer is capable of meeting the mixing requirements of the most popular consumables used in today's lab without the need to swap out head units. The truly innovative design makes maximum use of the space available, providing holding slots for no less than five different sizes of tube and a standard assay plate, while the traditional center cup can be used for quick hand-held mixing of a vast range of tubes and vessels.

- Variable speed rotary dial provides finite speed control during vortex mixing
- Stable elastomeric device feet, provide vibration dampening and prevent movement during operation
- Innovative head can be removed for cleaning purposes
- High impact ABS body shell

| Cat. No. | Description | Qty |
|----------|---------------------------------|------|
| 6100-800 | Vortexer, 110–240 VAC, 50/60 Hz | Each |


For pricing, see analytical-sales.com



Specifications:

- Speed Range: Variable 500-3000 rpm (shaking to vortexing)
- Orbit: 4.2 mm
- Control: 3 position switch—continuous, intermittent, manual, and off
- Weight: 8.5 Lb/3.8 Kg
- Dimensions: 173 W x 198 D x 198mm H (with Headpiece in place)
- SGS-US-C, CE, CB Scheme, RoHS 2, PSE, and WEEE compliant
- Two Year Warranty

Holding slots for no less than five different sizes of tubes and a standard assay plate

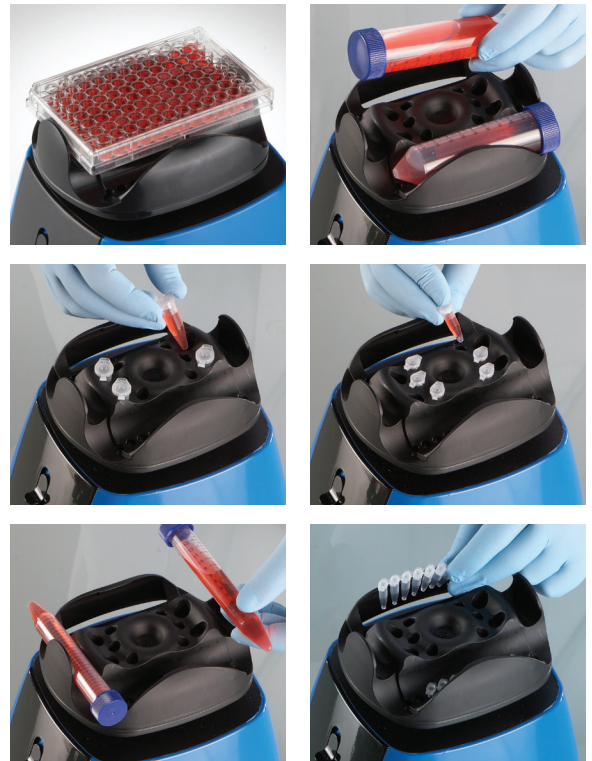

“All-in-one” head (which includes a center cup) can hold:

- Assay Plate x 1
- 50 mL conical tubes x 2
- 15 mL conical tubes x 2
- 1.5/2 mL microfuge tubes x 4
- 0.5 mL microfuge tubes x 6
- 0.2 mL microfuge tubes x 24

No Need to Change Heads!