



## **NEW! Gen II Column Heater Sleeve Controllers**

The Monocon and Multicon are the latest column heater controllers from Analytical. They are the next generation of the original MonoSleeve and MultiSleeve, and control Analytical's new Gen II HotSleeve+ column heaters. Significant new features include:

### **Touchscreen Interface**

- Provides simple and intuitive interaction with the device
- Provides greater utility to operators by allowing for gloved operation
- Tabbed, modern and familiar

### **Continuous Display of Critical Performance Data, Including:**

- Current temperature
- Set point temperature
- Elapsed time or time remaining until automated cutoff
- Activation status of heating unit
- Specifications and parameters of active heaters

### **High Power Operation**

- Support for heaters with a temperature range of 5°C above ambient to 130°C
- Faster heater to-temperature time allowing for decreased downtime
- Supports 250 W heaters allowing for increased performance and higher capacity operation
- Supports Peltier coolers

### **Screen Lock**

- Prevents accidental touches from changing settings

### **Stacking Enclosures**

- Preserves space by allowing multiple devices to share the same 16.5" x 11" vertical footprint

### **1000 Hour Maximum Timer Limit**

- Provides users with the time that they need to collect all necessary data

### **USB Connectivity**

- For PC connection

### **Modbus RS485**

- Allows the controller to be connected as part of a Modbus RS485 chain. This allows up to 31 controllers to be used simultaneously.
- A list of Modbus serial commands is available

Cat. No.	Description
Monocon	Monocon Gen II Single Column Heater Sleeve Controller
Multicon	Multicon Gen II Multiple Column Heater Sleeve Controller

## Monocon



Front




Back

## Multicon



Front

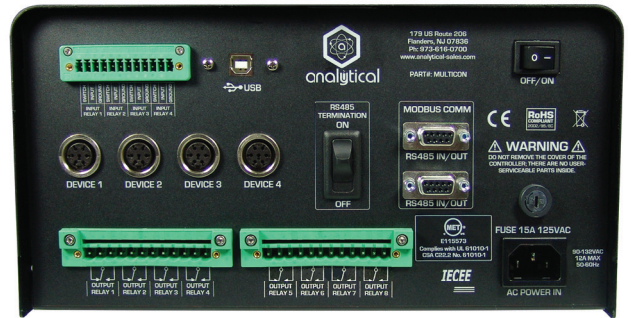

Back



## Gen II HotSleeve+ Self-Calibrating “Smart” Column Heaters

The Gen II HotSleeve+ is the second generation of the HotSleeve+ column heater, and the third iteration of the original Hot Sleeve. Use with either the new MULTICON™ Gen II multiple heater controller or the new MONOCON™ Gen II single heater controller, sold separately.

- Inner diameter: 0.65” – Improved internal geometry accommodates a superior fit for analytical columns
- Models designed to fit columns ranging from 5cm to 30cm long
- Temperature range: 5°C above ambient to 130°C
- Robust, high temperature VELCRO® fastener provides superior performance when compared to Gen I HotSleeve+ snap fasteners

Cat. No.	Description
HS2-5L	Gen II, HotSleeve+ Self-Calibrating “Smart” Heater for 5cm Columns
HS2-10L	Gen II, HotSleeve+ Self-Calibrating “Smart” Heater for 10cm Columns
HS2-15L	Gen II, HotSleeve+ Self-Calibrating “Smart” Heater for 15cm Columns
HS2-25L	Gen II, HotSleeve+ Self-Calibrating “Smart” Heater for 25cm Columns
HS2-30L	Gen II, HotSleeve+ Self-Calibrating “Smart” Heater for 30cm Columns